



GemLock™ Touchless Delivery System



© 2005 Zimmer Dental Inc. All rights reserved. 4903, Rev. 11/04.



Control. Confidence. Superior tactile response.

For more information about our Products, Professional Programs and Continuing Education, contact us:

1900 Aston Avenue
Carlsbad, CA 92008-7308, USA
In the U.S. 800 854 7019

To fax an order 888 225 2483
Outside the U.S. +1 760 929 4300
Australia +61 2 4297 4182 or 1800 241 916
Canada +1 800 265 0968 or 1 905 567 2073

France +33 (0)1 45 12 35 35
Germany +49 (0)761 4584 722/723
Israel +972 (0)3 6124242
Spain +34 93 846 05 43

To receive our eNews visit us at http://www.zimmerdental.com/news_eNewsLetterSignUp.aspx



www.zimmerdental.com



Precision you can trust. Performance that delivers.

In situations requiring extreme accuracy, there is no margin for error. As a gymnast grasps the rings, a firm, dependable grip is critical. Control, confidence, and a secure connection are essential for a perfect performance.

Take a close look at the *GemLock* Touchless Delivery System from Zimmer Dental. Utilizing advanced patent pending technology and a unique design, the *GemLock* system delivers a new level of performance and reliability with dynamic tactile response and touchless transfer during surgery. The result is a simpler, safer surgical experience with superior control and consistency.

Advanced ruby insert technology.

The *GemLock* delivery system utilizes spring-loaded synthetic ruby bearings to accurately and repeatedly secure the implant to the delivery tool.

Synthetic rubies possess remarkable strength and hardness, offering extreme abrasion resistance and low breakaway friction, resulting in expanded longevity and performance of your surgical tools. Manufactured in strict, controlled environments, these jeweled bearings also possess precision tolerances, and a smooth micro-finish for long-term consistency. Simply more dependable, the *GemLock* delivery system is the smarter alternative.



Feel the GemLock system difference.

GemLock delivery tools are designed using advanced technology and the highest quality materials to ensure consistent performance every time.


Enhanced reliability: *GemLock* touchless delivery tools offer a secure, reliable fit and superior control when interlocking components are mated.

Improved consistency: Durable synthetic ruby inserts improve longevity and ensure delivery tools remain consistent over time.

Internally engaging components: *GemLock* tools engage the Fixture Mount/Transfer internally for increased precision, making them ideal for use in tight interdental spaces.

Touchless transfer/delivery: Hands-free delivery from package to patient increases efficiency, decreases chair time, and minimizes the potential for contamination of the implant.

Proudly offering one of the most comprehensive dental implant product lines available, Zimmer Dental is a market leader in the development of world-class implantology products, practice-building strategies and educational programs focused on empowering clinicians and improving patients' lives.

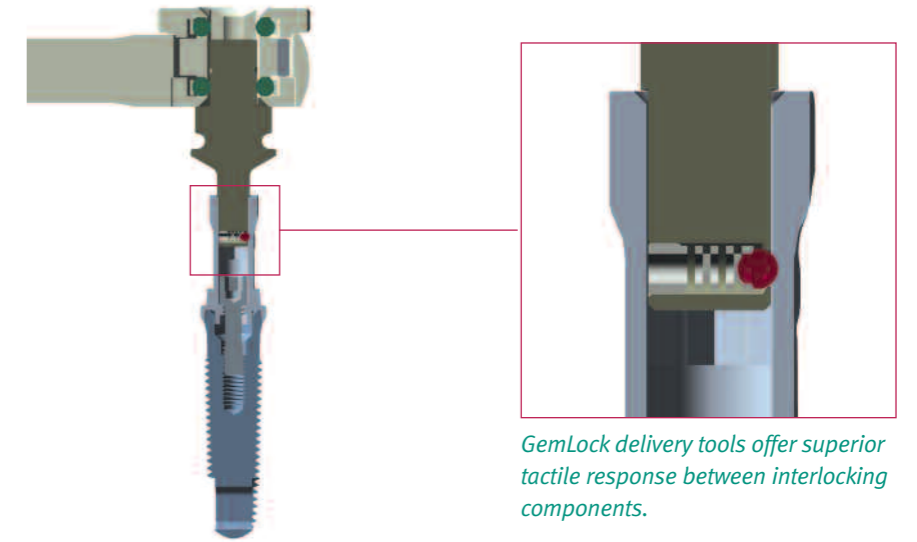

GemLock delivery tools offer superior tactile response between interlocking components.

GemLock Delivery Tools	Part Number
Retaining Square Ratchet	RSR
Retaining 2.5mm Hex Driver, Short	RH2.5
Retaining 2.5mm Hex Driver, Long	RHL2.5
Retaining 2.5mm Hex Driver, Latchlock	RHD2.5

***GemLock* Touchless Delivery System is available for Tapered Screw-Vent®, Screw-Vent®, AdVent®, and SwissPlus® implant systems. For more information or to order your set of delivery tools, call your local Zimmer Dental Sales Consultant or Customer Service at 800 854 7019.**