



## Locator® Overdenture Attachment



### Improving your patient's quality of life.

Retaining an overdenture with implants can greatly improve chewing function for the patient. With the *Locator* Overdenture Attachment from Zimmer Dental, patients are able to easily seat their overdenture, which will help prevent worn or destroyed attachments caused by improper seating. The *Locator* Overdenture Attachment is the ideal solution for overdenture restorations.

- **Low Profile.** A low vertical height makes the *Locator* Overdenture Attachment an effective option for treatment of a wider range of patients.
- **Easy for Patients.** Self-locating attachment gives patients the ability to easily position and place their overdentures.
- **Predictable.** Dual retention and tremendous durability ensure long-lasting stability and performance.<sup>1</sup>
- **Versatile.** A variety of options for retention and angulation correction let you tailor the case to fit the needs of the patient.


<sup>1</sup>Data on file.

The *Locator* Overdenture Attachment's unique combination of features offers outstanding advantages for you and your patients. The durable, low-profile design can be used in a wide range of overdenture cases and makes it easy for patients to seat their own denture correctly every time.

#### Low Vertical Height:

- Low-profile design provides a critical advantage in tight interocclusal spaces.
- Total abutment and attachment height is as low as 2.5 mm on a *Tapered Screw-Vent*® Implant.
- A variety of cuff heights allows you to accommodate varying tissue depth.

## The right solution for you and your patients.

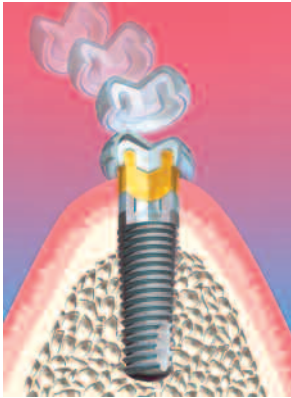

The self-aligning design of the *Locator* Attachment system helps guide the attachment into place on the abutment.



### Self-Aligning Design:

- Self-locating design helps guide the attachment into place on the abutment.
- Patients can easily position and place their overdenture every time.
- Reduces wear on the nylon liners due to improper seating.
- Longer-lasting liners mean fewer recall visits for the patient.<sup>1</sup>

### Incredible Durability:

- The *Locator* Overdenture Attachment retains on the inside and outside of the abutment for a stronger, longer-lasting fit.
- The unique dual retention design provides nearly twice the retentive surface area over internal retention attachments.
- The nylon male pivots inside the metal housing for a resilient connection.

### Versatility:

- Flexibility in angulation and retention with a variety of nylon inserts.
- Three levels of retention are available for cases up to 10° off axis per abutment: 1.5 pounds (blue), 3 pounds (pink) and 5 pounds (clear).
- Accommodates up to 40° of divergence for two implants (or 20° for one implant) with use of special nylon males with two levels of retention: 4 pounds (green) and 1.5 pounds (red).



*Locator* Overdenture Attachments are available for *Tapered Screw-Vent*, *AdVent*®, *SwissPlus*® and *Spline*® Implant Systems. For part numbers and pricing, please refer to the Zimmer Dental Price List for each product line.

To learn more about this product, please visit us online at [www.zimmerdental.com](http://www.zimmerdental.com) or to speak to a sales representative, call **800 854 7019**.

To receive our eNews visit us at [http://www.zimmerdental.com/news\\_eNewsletterSignUp.aspx](http://www.zimmerdental.com/news_eNewsletterSignUp.aspx)

For more information about our Products, Professional Programs and Continuing Education, contact us:

In the U.S. 800 854 7019  
To fax an order 888 225 2483  
Outside the U.S. +1 760 929 4300  
Australia +61 2 4297 4182 or 1800 241 916  
Canada +1 800 265 0968 or +1 905 567 2073

France +33 (0)1 45 12 35 35  
Germany +49 (0)761 4584 722/723  
Israel +972 (0)3 6124242  
Spain +34 93 846 05 43

Zimmer Dental  
1900 Aston Avenue  
Carlsbad, CA 92008-7308  
USA

[www.zimmerdental.com](http://www.zimmerdental.com)

