



Surgical Options Guide



Tapered Screw-Vent® Implant System –
Traditional and guided surgery instrumentation options.



Zimmer proudly introduces a comprehensive line of existing and future surgical instrumentation products that feature easily customizable configurations designed to meet your individual needs and surgical style.

The popular *Zimmer*[®] Instrument Kit System offers intuitive, color-coded drilling sequences for traditional surgery. For reduced chair time and clinician convenience, the soon-to-be introduced *Zimmer*[®] Drill Stop Kit is designed to take the drilling depth guesswork out of implant placement procedures.




Zimmer's flexible approach to guided surgery aims to allow placement of the flagship *Tapered Screw-Vent* family of implants with model- and software-based guides. Zimmer instrumentation is designed to enable implant placement with widely used treatment planning software systems.

Committed to helping you maximize your time and streamline your surgical procedures, Zimmer presents an expanded portfolio of surgical solutions with the ideal combination of flexibility, simplicity and convenience.

Zimmer Dental surgical product options

traditional surgery

traditional surgery

Surgical need	Traditional	Traditional with depth control
<p>Required instrumentation</p>	<p>Zimmer Instrument Kit System</p> 	<p>Zimmer Instrument Kit System*</p>  <p>Zimmer Drill Stop Kit*</p> 
<p>Position and angulation control</p>	<p>Visual assessment</p>	<p>Visual assessment</p>
<p>Depth control</p>	<p>Drill depth markings</p>	<p>✓ Yes</p>
<p>Requires treatment planning software</p>	<p>No</p>	<p>No</p>
<p>Approximate availability</p>	<p>Available now</p>	<p>Early 2010</p>

* Requires updated *Driva™* Drills (marked with axial stripes)

Zimmer Dental surgical product options

guided surgery

Model-based surgical guide

Software-based surgical guide

Software-based surgical guide

Zimmer Instrument Kit System*




Zimmer Drill Module with *Dríva* EG Drills



Zimmer Instrument Kit System*



Zimmer Drill Module with *Dríva* EG Drills

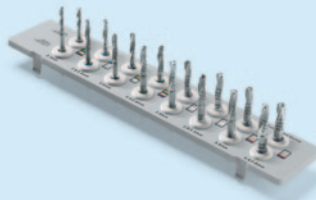

Zimmer® Tube Adapter Kit



Stay tuned for more advancements in the Zimmer Guided Surgery Portfolio.



For drilling only

For drilling only

No

✓ Yes

No

✓ Yes

Early 2010

Early 2010

2011

guided surgery

* Requires updated *Dríva* Drills (marked with axial stripes)

Zimmer Dental surgical product description


traditional & guided surgery



Part No: TZKIT



Part No: DSKIT

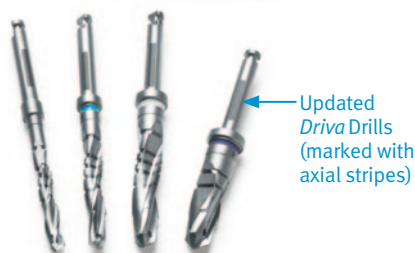


Updated
Driva Drills
(marked with
axial stripes)



Part No: GSMOD



Part No: TADKIT

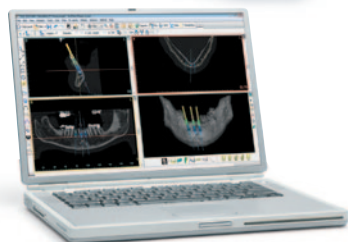

Zimmer Instrument Kit System

From complete set-ups that include all instruments, to standalone instrument kits and a unique Staging Block, the *Zimmer Instrument Kit System* is conveniently adaptable to your individual needs. Intuitive instrument organization and color-coding make the surgical sequence easy to learn and follow.

Zimmer Drill Stop Kit

The new Drill Stop Kit includes a set of titanium, reusable drill stops designed to limit drilling depth from bone level during osteotomy preparation. Featuring a convenient “pick and go” stop application mechanism, this cost-efficient kit is designed to save chair time and increase clinician convenience. Drill Stops are only intended for use with [updated Driva™ Drills \(marked with axial stripes\)](#).

Updated Driva Drills

Designed with a reduced hub, flange and two axial stripes 180° apart for easier identification, these updated drills are required for the use of the *Zimmer Drill Stop Kit* and/or guided surgery instrumentation. These drills will replace the existing line of *Driva Drills*, starting in December 2009 (note: part numbers will remain the same).

Driva EG Drills

Additional *Driva EG Drills* are available with longer shafts, in two new lengths.

Zimmer Drill Module

This kit insert includes sixteen *Driva EG Drills* and can be snapped into your *Tapered Screw-Vent Surgical Kit*.

Zimmer Tube Adapter Kit

Designed to fit in the tubes located inside model- and software-based surgical guides, these surgical instruments orient drills and provide positional and angular control.

Zimmer Guided Surgery

Stay tuned for more advancements in the *Zimmer Guided Surgery Portfolio*.

Treatment Planning Software

Widely used treatment planning software and guided surgery systems.

Please contact a Zimmer Dental representative for more information.

For more information about our Products, Regenerative Materials and Educational Opportunities, contact us:

1900 Aston Avenue
Carlsbad, CA 92008-7308
USA

In the U.S. 1 (800) 854-7019
To fax an order 1 (888) 225-2483
Outside the U.S. +1 (760) 929-4300
Australia +61 (0)2 9950 5444 or 1 (800) 241-916
Canada +1 (905) 567-2073 or 1 (800) 265-0968

France +33 (0)1 45 12 35 35
Germany +49 (0)761 1 56 47 0
Israel +972 (0)3 6124242
Italy +39 043 855 5573
Spain +34 93 846 05 43

To receive our eNews visit us at http://www.zimmerdental.com/news_eNewsLetterSignUp.aspx



www.zimmerdental.com